

# KL 1500 HAL

A benchmark for colors and functionality



KL 1500 HAL

The KL 1500 HAL is a 150 W halogen cold light source, developed and designed for industrial and life science applications in stereo microscopy and macroscopy.

Close to 100 CRI halogen spectrum for excellent performance at rendering colors.

The optical system is optimized for coupling the light of the EFR-type halogen lamp into SCHOTT KL 1500 light guides.

General data			
Part-No.			150 700
Dimensions (W x D x H)	(mm) / (inch)		200 x 286 x 147/7.9" x 11.3" x 5.8"
Weight	(kg)		approx. 4.2
Cooling			ultra low-noise fan
Electrical data			
Operating voltage, frequency			100 – 240 V ~, 50 – 60 Hz
Power consumption	(VA)		180
Protection class			I
Overvoltage category			II
Lamp type			Halogen reflector lamp type EFR
Lamp voltage rating	(V)		15
Lamp power rating	(W)		150
Average Lamp	ECO	(h)	1500
Service life	BOOST	(h)	150
Illumination data			
Light flux (Ø 9.0 mm / 1000 mm P / N 155104)	(lm)		600 <sup>1</sup>
Light control			Electrical
Active light guide diameter	(mm) / (inch)		Max. 9/max. 0.35"
Certification marks			
Marking			CE
Approvals			CS US
EMV-class			B

Note: country specific power cord must be ordered separately!

<sup>1</sup> at the output of SCHOTT fiber optic light guide, 1-arm, flexible, Ø 9 mm active, length 1000 mm





## Features:

- Internal wide range power supply (100 – 240 V, 50 – 60 Hz)
- High frequency illumination (20 kHz, digital camera proof)
- Continuous dimming from 0 to 100 %
- Swing-in optics for homogenization of light output
- High quality housing with integrated handle
- Ultra low noise fan
- Filter slider with two positions (filter can remain in the light source even if it is not used)
- Lamphousing opens to the side (allows change of bulb even if stacked)
- Easy to use fixture mechanism with clamping screw

Power cords	Part-No.
IEC cord with European plug	400 051
IEC cord with US plug	400 052
IEC cord with UK plug	400 053
IEC cord with Swiss plug	400 054



KL 1500 HAL

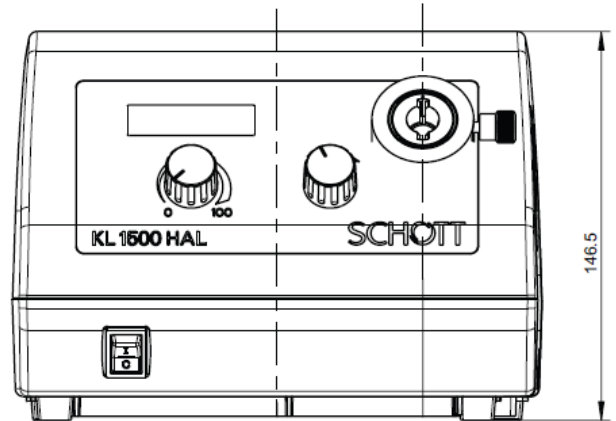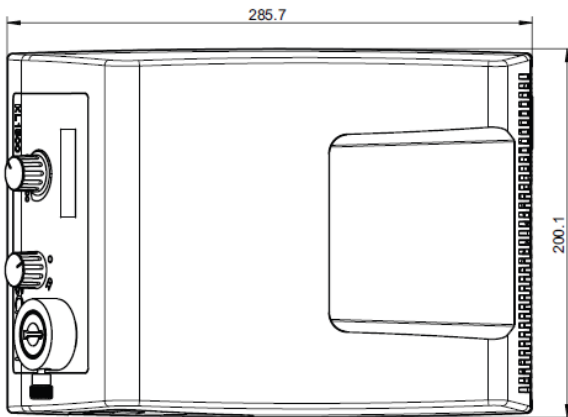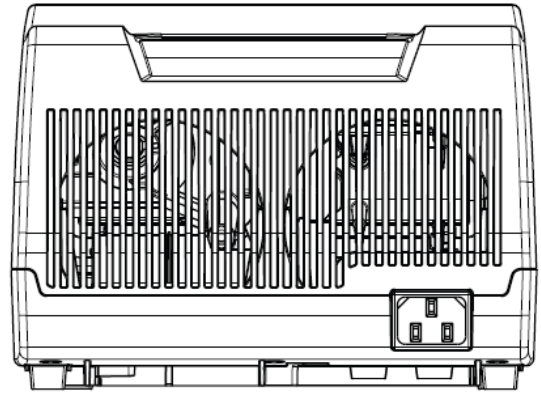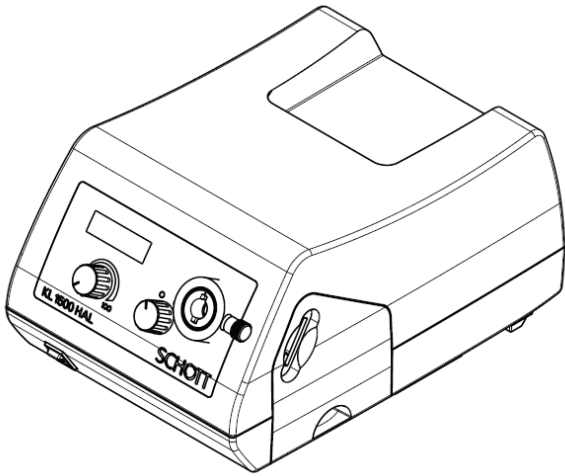

All specifications are subject to change without prior notice. This datasheet or any extracts thereof may only be used in other publications with express permission of SCHOTT. © SCHOTT North America, Inc.

SCHOTT North America, Inc.  
**Lighting and Imaging**  
 122 Charlton Street  
 Southbridge, MA 01550 USA  
 phone: +1-508-765-3211  
 fax: +1-508-764-6273  
 lightingimaging@us.schott.com  
 www.us.schott.com/lightingimaging

